



# CROSS POINT

## Machines

Laser

RF

Facial Machines

Body Machines

Carboxy Therapy (Co2)

Meso Gun (with-Without) Needles

## Products

Stem Cell Therapy

Mesotherapy

Consumables

Retail

# CO2 FRACTIONAL

Laser Type:	10.6μm
Output Power:	40W
Focal length of working head:	F=50mm 75mm 100mm
Indication Beam:	Red Diode Laser(650nm, ≤5mW)
Transmission System:	7-Articulation Joint Spring Arm
Smoke Exhausting System:	Built-in air blow
Cooling System:	Air cooling
Environmental temperature:	5° C~40° C
Relative Humidity:	≤80%
Power Supply:	~220V ± 22V, 50Hz ± 1Hz
Dimension:	52*68*132cm
Weight:	35kg

CO2 laser treatment uses a pulsed emission laser to form a burning zone composed of laser action points and spaces in the epidermis. Each laser action point consists of a single or several high-energy laser pulses that can penetrate directly into the dermis layer. Instantly vaporizes the tissue of the lesion, stimulates the proliferation of collagen, and initiates a series of skin reactions such as tissue repair and collagen rearrangement. The patient can return to normal life after four to five days. The treatment itself is safer and can treat any part of the body.

The product has novel design, stable work, short irradiation time and long-lasting effect; easy to master, good clinical effect and high safety factor.



# Q-SWITCH LASER


Energy Level:	Single mode: Max 400mj (532nm), Max 800mj (1064nm) Double Mode: Max 800mj(532nm), Max 1600mj (1064nm)
Wavelength:	1064nm, 532nm
Pulse Duration:	750ps
Repetition frequency:	1-10Hz
Spot size:	2-10MM
Lamp Quantity:	1 lamp 1 rod / 1 lamp 2 rod
Dimension:	90×45×170cm
Weight:	105kg



## Laser source: Q-Switch ND: YAG solid-state laser

Wave length:	1064nm, 532nm
Output Energy:	Single pulse: 1064nm: 100-800mj; 532nm:100-400mj
Energy density (max.):	Double pulse: 1064nm: 100-1600mj; 532nm:100-800mj
Pulse duration:	1064nm: 15.0J/cm2; 532nm:8j/cm2
Repetition Rate:	750ps
Spot Size:	1-10hz adjustable
Light guide beam:	φ2-10mm adjustable
Cooling System:	650nm red laser
Laser delivery:	Closed Cycle Water to Air Heat Exchanger
Controlling system:	Korea 7 jointed articulated arm
Power requirement:	10.4inch TFT color touch screen
Machine dimensions:	220V,10A,50/60Hz; 110v,10A, 50/60Hz
Net weight:	680mm×400mm×820mm
Note:	75kg
Technical data changed will not inform specially	


# Electrolysis

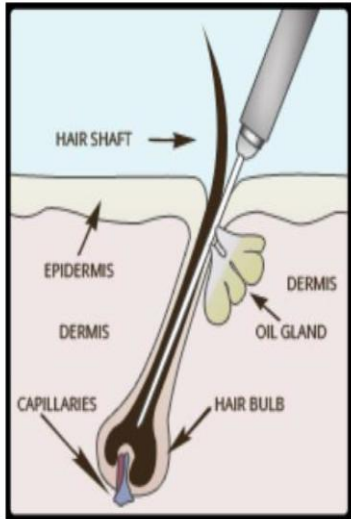

Electrolysis is a method of removing individual hairs from the face or body areas of the body can be treated with electrolysis, and there are general permanent side effects.



Install the Needle

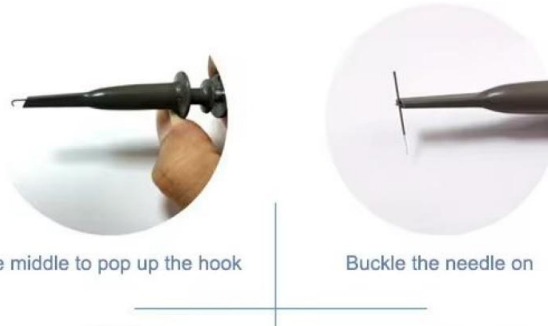

Press the middle to pop up the hook

Buckle the needle on



Specification:

Input power :500 W

Output power :0-100%  $\pm$  10 % adjustable

Output Frequency :13.26MHz

Output Voltage :AC 30-150V

Needle diameter :0.01mm

Output mode :Continuous or pulse for option

Display :8.4 inch

Function :Vascular removal, Hair removal, white hair removal


Operation method:

- ◆ Plug the accessory foot switch into the output jack on the front of the instrument.
- ◆ Connect the net wire pen to the host.
- ◆ Disinfect the skin and needle tip, then turn on the power, select the output gear and select the appropriate gear for operation.
- ◆ Aim the needle tip at the blood supply point, or the main trunk, to reduce pressure, and do not pierce the skin. (It is normal that there will be liquefaction precipitation) until the blood vessels disappear, the interval between every two points should not exceed 1cm, and the time should not be too long at one point.
- ◆ After using, turn off the power and unplug the power plug
- ◆ Use with caution for people with scars. Large blood vessels or dense blood vessels often cannot disappear completely at one time. If a second treatment is required, it must be performed 15 days apart before operation. It is not advisable to get wet within three days to prevent infection.

# HIFU 10 LINES (30K)

## • Face Lifting (3 Cartridges)

4D HIFU works for face lifting, skin tightening, wrinkles, pores, V-line forming, jowl line, cheek, skin tone improvement, and laxity improvement.

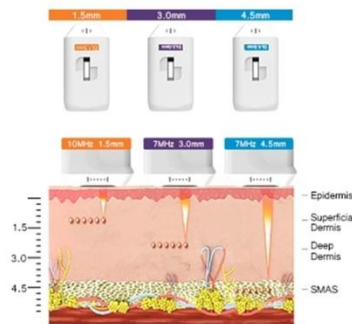

## • Body Slimming (5 Cartridges)

Whole body contouring of chest, arm, abdomen, flank, thigh, lower leg, and back.



### Type

Frequency

Power

Treatment Depth

Treatment Length

Width

Energy

Pitch

Clearance

Power Source

Machine Size & Net weight

Package Size & Gross weight

Trolley Size & Net weight

Trolley package Size & Gross weight

### High Intensity Focus Ultrasound

4MHz

200W

DL1.5mm DL3.0mm DL4.5mm (Standard)

DL6.0mm DL8.0mm DL10.0mm DL13.0mm

DL16.0mm (Optional)

5 -- 25mm

0 -- 10mm

0.2 -- 2.0j

1.0mm-5.0mm

1.0mm-5.0mm

AC 110V~220V~50Hz/60Hz

400X350X335mm & 8.5Kg

525X445X465mm & 12Kg

540X480X790mm & 10Kg

680X680X330mm & 15Kg

# CARBOXY THERAPY (CO2)

Screen size	10 inch touch screen
Programming modes	Automatic system and Manual system to optional
Function	Acne, lose weight, <u>Cellulitis</u> , stretch marks, hair growth, <u>dark undereye circle</u> , etc.
Temperature	Automatically heats when the temperature is below 23 degrees.  Stop heating when the temperature is above 30 degrees. (output gas measurement)
Gas in Pressure	4.00-10.00bar
Gas in pressure	From 10-150ml/min
<u>Biosafety filter</u>	0.22 $\mu$ Sterile filters
Recommended needles	30g*1/2 or 27g*1/2
Machine size	420*330*250mm
Packing size	565*475*425mm
Net Weight	10 Kg
Voltage	110/220V
Delivery parts free	High pressure regulator 1 set  Input gas connecting pipe 1 piece  Output gas connecting pipe 1 piece  Foot pedal 1 set  Sterile filter two pieces  Aseptic needles 4 pieces  Fuse two pieces  Aseptic extension tube 2 piece





# EMS SCULPT & SLIMMING

- Magnetic wave
  - Voltage
  - Out put power
  - Frequency
  - pulse width
  - Mode
  - screen
  - Treatment handle
  - machine size
  - package size
  - Net weight
  - Gross weight
- 0-7tesla
  - 220V
  - 4300W
  - F1:1-10Hz F2:1-100Hz
  - 300us
  - model-I (Sport mode)
  - model-II (Gentle mode)
  - model-III (Professional mode)
  - 10.4 inch
  - 4 handles
  - 1500mm\*360mm\*520mm
  - 1150mm\*540mm\*840mm
  - 75kg
  - 105kg





# 360 CRYO LIPOLYSIS

- Input power: 1200W
- Frozen Head Pressure: 0-60Kpa
- Frozen Head Screen: -12°C -1°C
- Treatment Handle's Type: 5 applicator in different size
- Machine Screen: 15.6 inch Collapse screen, 30°-73° adjustable.
- Cooling System: Semiconductor cooling + water cooling +air cooling
- Dimension: 41\*58\*154cm

After the adipose tissue of the target care site is sucked into the cooling probe, the cooling probe releases a 360-degree surrounding cooling energy to act on the fat cells without damaging any surrounding tissue. Cryocool's steady and powerful cooling energy triggers a natural cell death process called "apoptosis", after which apoptotic fat cells are excreted through the body's natural metabolic process. Cryocool targets fat cells for permanent removal from the body, resulting in significant fat loss in the treated area.

## ADVANTAGES



Medical silicone, safe and comfortable



Medical filter cotton, effectively filter debris



# FACIAL MULTI ARMS

## ► One Machine Multi-functions


Ⓐ

**Reskin Handle**  
Patented Gas-liquid  
Infiltration Technology ,  
Nano-nutrient absorption

With 450m/s speed and 50-850Kpa pressure  
jet technology, deep hydration and  
long-lasting effects, and  
stimulate collagen




Traditional Needle Injection

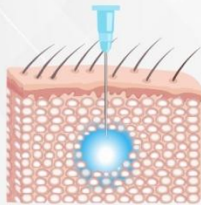


(20 days)

Needle Free Injection

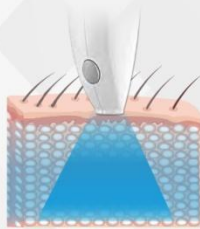


(3-7day)

Wrinkle production →  
As age increasing,  
collagen and elastin will lost  
and skin get wrinkles

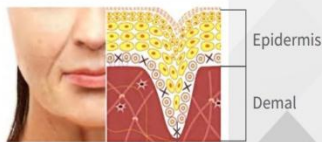

Protein regeneration →  
With needleless injection,  
more even distribution

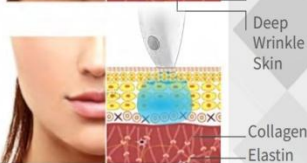

## THE HANDLE



Reskin Handpieces



Hydro Dermabrasion



Skin Scrubber



RF lifting



Ultrasonic



Cold Hammer

## SPECIFICATION

Power: 500W

Adjustable Mode: CON&INT

High Pressure: 5KG/CM<sup>2</sup>-10KG/CM<sup>2</sup>; 850Kpa

Penetrate Depth: 0.2MM-0.3MM

Dose capacity: 1UL-24UL

Cold Hammer Temp: 0-5°C

Wet Standing Vacuum Flow: 30L

RF Output : 5MHZ

Ultrasonic: 1M



# MESO TWIN GUN (WITH-WITHOUT) NEEDLES

- Acne problems
- Congestion
- Blemishes
- Dull & sagging skin
- Dehydration
- Fine lines
- Sensitive skin



## » Expected Effect

1. Improve dull skin to clear, flawless and transparent skin
2. Smooth skin texture
3. Tighten pores by removing wastes in pores
4. Control sebum secretions
5. Absorb cosmetics quickly and take make-up well



# U.S. CAVI MULTI ARMS

## Technical specifications



Handle type	Frequency	Treatment energy level	Treatment mode	Tip size
Focus ultrasound	230Khz	1-30 degrees	Pulse mode	D42mm
Ultrasound cavitation	40khz	1-10 degrees	Continuous mode	D50mm
Bipolar RF	2mhz	1-30 degrees	Continuous mode	D8mm, Distance:15mm

Input Power	10-400w
Focus Depth	12mm
Cooling head temperature	-5℃—5℃
Treatment life span of hand piece	500,000 shots
Machine dimensions	550X360X1175mm
Package dimensions	705*515*1130mm
Machine weight	38kg
Package weight	63kg
Power requirements	110V-240V,10A single phase
Environmental requirements	+15℃—30℃, 80% humidity

4





# MIRCONEEEDLE RF MULTI ARMS

## MICRO NEEDLE RF Applicator

- MICRO NEEDLE RF (RF) applicator for procedures requiring ablation and resurfacing of the skin. Add vacuum  
- suction combined probe for better skin contact, especially around eyes and neck areas.



M-18

### Stepping Motor Type

- Stepping Motor Type, different to existing solenoid type, the needle inserts to the skin smoothly with no shock, and causes no bleeding and no pain after the procedure.



M-49

### Gold Plating

- Needle is durable and also have high biocompatibility by applying Gold Plating. Patient with Metal allergy could also use it with not concerning Contact Dermatitis.



M-25

### Needle Depth Control

- 0.4-4.0mm (0.2mm Step)
- Operate epidermis and dermis layer by controlling the needle depth in unit of 0.2mm.



M-81

### Safety Needle System

- Sterilized disposable needle tip
- operator can easily notice the applying RF energy from LED light.

### Needle Thickness: 0.2mm

- The needle structure is easy to penetrate the skin with minimum resistance.



Power	max 50 watts
Handle	Fractional RF handle Micro Needle RF handle
Frequency	2MHz
Energy	1-10 Level
Pulse Duration	10-900ms
Needle Depths	0.4-4mm adjustable
Pressure Level	1-4
Working mode	SINGLE/AUTO
Mode	MONO & BIPOLAR
Cooling	1-3 Level (Max -10°C)
Tips	Insulated & Non-insulated
Invasive Tips	M-18, M-25, M-49, M-81
Non invasive Tips	F-3*3 , F-5*5, F-9*9
Dimension	46*37*125CM
Net weight	
Power	10% , 10A , 50HZ / AC110V±10%10A , 60HZ

# PLASMA PEN

Items	Ultra V Plasma	Others
<b>Tip</b>	* Non-toxic raw material (SUS 316 – does not contain lead)	* Some part of risky materials (Copper – contains lead Iron – tetanus (lockjaw) possibility hazard)
<b>Power</b>	* Sparkling technique ✓ Does not touch the skin, minimal pain ✓ Shortens treatment time, due to higher power	* Scratching technique ➢ Severe pain ➢ Severe damages depending on skill
<b>Adaptor</b>	* 2.5m length ✓ Comfortable to use ✓ Easy to treat even hard to reach places	* 1-1.5m length ➢ Not comfortable to use due to the short cord
<b>Packing</b>	* Fine pouch ✓ Comfortable to store and carry (no need to purchase a case separately) ✓ Luxurious material – even though there are many third party cases available, none can compare to the material, feel and design of this pouch ✓ Quality and appearance oriented clients choose Ultra V Plasma	* Carton box ➢ Uncomfortable to carry ➢ Appearance is not luxurious ➢ Case needs to be purchased separately
<b>Safety</b>	* Stable power * Dual safety device	* Break or heat after 10 min of use * Emit radiofrequency right after turn-on * No safety function
<b>Certificates</b>	* KC * CE * FCC	* NO



Part	Power	tip	remarked
Eyelid lift	Power 2-3	T1	Usual method
Double eyelid	Power 2-3	T1	Usual method
Fine wrinkle	Power 1-2	T1, T2	Depends on the condition
Deep wrinkle	Power 3-4	T2/T3	
Pigment spot	Power 2-3	T1/T2/T3	
Scar & stretch marks	Power 2-3	T1/T2/T3	

- Frequency : 2-6 treatments in one month.
















# CHAIR



# VELA FAMILY

PRODUCTS			
CLASSIFICATION	PIC.	PRODUCT NAME	SPEC
VELA		VELA CONTOUR V COTNOURING SERUM FOR JAW AND NECKLINE	30ml spoid dropper
		VELA CONTOUR 4D V REFILL PACK	10ML X 55SHEETS
		VELA CONTOUR V FIRMING CREAM	30g
		VELA CONTOUR 4D V LIFT BAND	VELCRO TYPE BAND
		VELA CONTOUR ULTRA V LIFTING BAND	VELCRO TYPE BAND
		VELA CONTOUR 4D FIRM MASK (25STEP)	24gX5pcs , 7g FIRMING CREAM
		VELA CONTOUR DUAL LIFT MASK	10.5gX5pcs
		VELA SHAPE RESCULPTING GEL	150ml Cream type
		VELA SHAPE REAFFIRMING GEL	150ml Cream type
		VELA SHAPE ACTIV CREAM	150ml Cream type
PHOSTILYN		INJECTION METHOD	10X10ml Vial
Drug		ANESTHESIA CREAM (LIDOCAIN 9.6%) - DRUG	500G
DERMATO COSMETOLOGY		K OXIDERM OP CREAM (VIT. K1 OXIDE 2%)	15ml tube
			30ml tube

CLASSIFICATION	PIC.	PRODUCT NAME	SPEC
LITHEMIS		RESYNC SERUM	20ml dropper
		5C INTENSIVE CREAM	150ml tube
ACLEON		ACLEON SEBODERM SOLUTION	15ml spoid dropper
		ACLEON SEBODERM SOLUTION	55ml spoid dropper
		ACLEON SEBO X TRACTOR	200ml
		ACLEON SEBO X TRACTOR	1000ml
		ACLEON ADVANCED SPOT CONCENTRATE	30g
		ACLEON ADVANCED SPOT CONCENTRATE	30g
ACITVE PPC		ACLEON SEBODERM CLARIFYING TONING LOTION	200ml
		ACLEON SEBODERM PURIFYING CLEANSER	200ml
		ACTIVE PPC CREAM	100g , tube
			500g , pump bottle
		ACTIVE PPC SOLUTION	100ml pump type



DUBAI, AL GARHOUD, RED AVENUE, OFFICE NO.1-18. | PHONE:  
+971 52 450 8255